


NEW!

# CliP® Piston MVP

*Titanium Piston for malleovestibulopexy with a unique ball joint for intraoperative adjustments, designed by Hausler.*

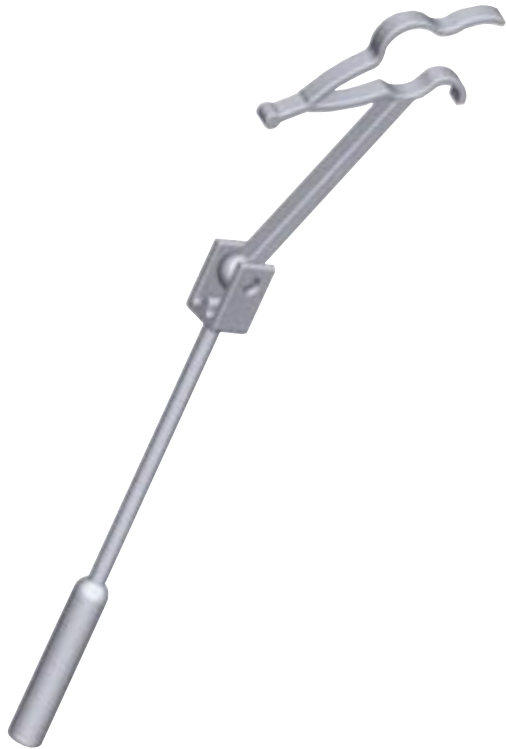

Stapes prostheses are used to replace fixated stapes. Normally the prosthesis is attached to the long process of the incus, but with an extra long shaft, the prosthesis can be coupled to the malleus handle. The stapes prosthesis then creates a direct connection between the malleus handle and the inner ear.

In collaboration with the Cantonal University Hospital of Berne, Switzerland and on basis of the CliP® Piston àWengen, a unique implant was developed with favourable intraoperative handling.

Until now, a normal loop piston was used when creating a bridge between the malleus handle and the inner ear. To get a satisfying coupling between the prosthesis and the malleus handle, the loop piston had to be bent in an almost right angle and fixed to the malleus handle by crimping.

With this unique Hausler designed stapes prosthesis, both the complicated bending and the critical crimping are eliminated. The bendable shaft with the new ball joint gives the CliP® Piston MVP an optimal intraoperative adjustment and the CliP® mechanism secures a solid connection to the malleus handle without the need of crimping.


## Benefits:

- New ball joint for optimal intraoperative adjustments.
- Proven CliP® design for a secure, crimpfree connection to the malleus handle
- Titanium implant with excellent biocompatibility
- Light weight for an optimal sound conduction
- MRI compatible



Easy bendable shaft due to ball joint.

Material: Titanium  
(ASTM F67 Medical Grade)  
Diameter: 0.4 / 0.6 mm

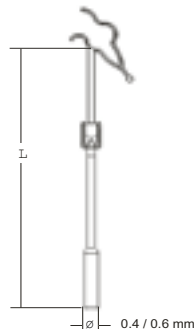

Length: (L) Ø 0.4 mm	Reference number	Length: (L) Ø 0.6 mm	Reference number
5.75 mm	1006 711	5.75 mm	1006 761
6.00 mm	1006 712	6.00 mm	1006 762
6.25 mm	1006 713	6.25 mm	1006 763

Hear the quality!™

© Atos Medical AB, 2004/06, REF 7834. Hear the Quality is a pending trademark owned by Atos Medical AB.

**ATOS**  
MEDICAL

CE 0124



Manufactured by:  
Heinz Kurz GmbH Medizintechnik, www.kurzmed.com



**International Marketing Partner:**  
Atos Medical AB  
P.O. Box 183, SE-242 22 Hörby, Sweden  
Tel +46 (0)415 198 00, Fax +46 (0)415 198 98  
info@atosmedical.com, www.atosmedical.com

**UK Distribution:**  
Platon Medical Ltd  
P.O. Box 2568, Eastbourne, East Sussex, BN21 3HZ, UK  
Tel 01323 431 930 • Fax: 01323 732 606  
E-mail: platon@mistral.co.uk